

UltraClean 3

Washer Disinfector Dryer



UltraClean 3



BenchTop



Under Bench

Technical Specification

Product Name: UltraClean 3 Washer Disinfector Dryer

Part No.	88002011	8812011
Name	Standard Bench Top	Standard Under Bench
Dimensions		
Width	560mm	560mm
Height	555mm	800mm
Depth	575mm	575mm
Max net weight	43kg	53kg
Gross weight	46kg	56kg
Chamber width	400mm	400mm
Chamber depth	450mm	450mm
Chamber height	300mm	300mm
Load capacity		
Chamber volume	60 litre	60 litre
Basket volume	50 litre	50 litre
Process monitoring	✓	✓
Fault alert	✓	✓
Warranty		
Parts (unlimited usage, excludes consumables)	1 year	1 year

Part No.	88002011	8812011
Name	Standard Bench Top	Standard Under Bench
Connections		
Drain connection	25mm dia	25mm dia
Voltage	230 VAC	230 VAC
Wattage	2750 W	2950 W
Frequency (hz)	50 Hz	50 Hz
Current	12A	12A
Power pump	150 W	150 W
Hydraulic connection cold water	¾" male	¾" male
min/max pressure	2.5 Bar min 5 Bar max	2.5 Bar min 5 Bar max
min/max temperature	+5 to +15 degrees C	+5 to +15 degrees C
Hydraulic connection warm water	¾" male	¾" male
min/max pressure	2.5 Bar min 5 Bar max	2.5 Bar min 5 Bar max
min/max temperature	+5 to +15 degrees C	+5 to +15 degrees C
maximum capacity	8 kg	8 kg
Noise	52 dB(A)	52 dB(A)
Max room temperature permitted	40 degrees C	40 degrees C

UltraClean 3

Bench Top & Under Bench Washer Disinfector Dryer

Proven performance and reliability made the UltraClean II range a popular choice for dental practices.

Now the new UltraClean 3 range brings the same high quality with the additional advantages of faster cycle times, a glass viewing panel and interior chamber illumination – all at a very competitive price.






More features and benefits:

- Incorporates a new, 35 minute 'Quick Cycle' for faster processing of instruments
- Viewing window and illuminated chamber enables the operator to make a simple check to ensure the machine is operating correctly
- Seamless hands free process from the dentist, through cleaning to sterilization in the autoclave, eliminating sharps injuries and instrument damage
- Fully compliant to EN 15883 and HTM2030
- High load capacity - despite its compact size, UltraClean 3 has one of the highest capacities on the market - for cost effective cleaning
- Forced hot air drying removes all moisture on instruments making them sterilization-ready for the autoclave
- Self-dosing detergent delivers consistent results - reduces waste and costs less per cycle
- Unique air gap prevents backflow pollution to ensure the safe disposal of contaminated fluid as required by EN 1717
- Can be supplied with optional printer or data logger for easy, recordable data capture
- Thermo-disinfects at 93°C for guaranteed cleaning results and features a built-in water softener



Furniture and Accessories

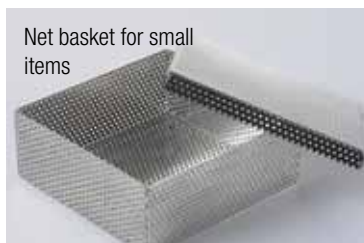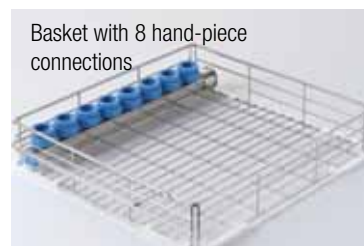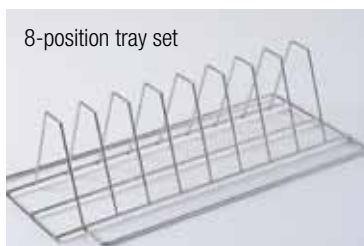

Prestige Medical - integrated decontamination solutions

For over 20 years, Prestige Medical has been respected as a class leader in the development and manufacture of sterilization and disinfection equipment. With infection control standards becoming increasingly stringent, a new breed of equipment is now required that is capable of meeting those challenges. The UltraClean 3 washer disinfectant helps you surpass those standards.

standards and guidelines

Designed and manufactured to comply with the following standards and guidelines:

- 93/42/EEC Medical Device Directive
- EN15883-1
- HTM 2030



Contact us for advice on infection control



Prestige Medical Limited
East House, Duttons Way,
Shadsworth Business Park,
Blackburn BB1 2QR

Tel : +44 (0) 1254 682 622
Fax : +44 (0) 1254 682 606

www.prestigemedical.co.uk
sales@prestigemedical.co.uk

@PrestigeMedica1

Registered in England
Reg No. 2826793
Vat No. 633707835

Registered Office:
Barlow Road,
Aldermans Green Industrial Estate,
Coventry CV2 2LD

Prestige Medical is a registered trademark



Due to a continuous process of development, Prestige Medical reserves the right to amend the specification of products and services without notice.