

# UltraClean 3 Washer Disinfector Dryer



Bench Top



## UltraClean 3

Bench Top & Under Bench Washer Disinfectant Dryer

Designed to meet 'best practice' requirements

Proven performance and reliability made the UltraClean II range a popular choice for dental practices needing to upgrade to meet best practice standards as defined by the HTM 01-05 guidelines.

Now the new UltraClean 3 range brings the same high quality with the additional advantages of faster cycle times, a glass viewing panel and interior chamber illumination – all at a very competitive price.

More Features and benefits...

- Incorporates a new, 35 minute 'Quick Cycle' for faster processing of instruments
- Viewing window and illuminated chamber enables the operator to make a simple check to ensure the machine is operating correctly
- Fully compliant to EN 15883 and HTM2030 - the recommended 'best practice' standard – to give you real peace of mind
- Seamless hands free process from the dentist, through cleaning to sterilization in the autoclave, eliminating sharps injuries and instrument damage
- High load capacity - despite its compact size, UltraClean 3 has one of the highest capacities on the market - for cost effective cleaning
- Forced hot air drying removes all moisture on instruments making them sterilization-ready for the autoclave
- Self-dosing detergent delivers consistent results - reduces waste and costs less per cycle
- Unique air gap prevents backflow pollution to ensure the safe disposal of contaminated fluid as required by BS EN 1717:2001\*
- Can be supplied with optional data logger or printer for easy, recordable data capture – essential for compliance
- Machines are available to meet specific Scottish NHS and England & Wales NHS requirements – please contact us for more information

\* Protection against pollution of potable water in water installations and general requirements of devices to prevent pollution by backflow



# UltraClean 3



BenchTop

Under Bench

## Technical Specification

Product Name: UltraClean 3 Washer Disinfectator Dryer

Part No.	88002110	88012110
Name	Standard Bench Top	Standard Under Bench
<b>Dimensions</b>		
Width	560mm	560mm
Height	555mm	800mm
Depth	575mm	575mm
Max net weight	43kg	53kg
Gross weight	46kg	56kg
Chamber width	400mm	400mm
Chamber depth	450mm	450mm
Chamber height	300mm	300mm
Load capacity		
Chamber volume	60 litre	60 litre
Basket volume	50 litre	50 litre
Process monitoring	✓	✓
Fault alert	✓	✓
<b>Warranty</b>		
Parts (unlimited usage, excludes consumables)	1 year	1 year


Part No.	88002110	88012110
Name	Standard Bench Top	Standard Under Bench
<b>Connections</b>		
Drain connection	25mm dia	25mm dia
Voltage	230 VAC	230 VAC
Wattage	2750 W	2950 W
Frequency (hz)	50 Hz	50 Hz
Current	12A	12A
Power pump	150 W	150 W
Hydraulic connection cold water	¾" male	¾" male
min/max pressure	2.5 Bar min 5 Bar max	2.5 Bar min 5 Bar max
min/max temperature	+5 to +15 degrees C	+5 to +15 degrees C
Hydraulic connection warm water	¾" male	¾" male
min/max pressure	2.5 Bar min 5 Bar max	2.5 Bar min 5 Bar max
min/max temperature	+5 to +15 degrees C	+5 to +15 degrees C
maximum capacity	8 kg	8 kg
Noise	52 dB(A)	52 dB(A)
Max room temperature permitted	40 degrees C	40 degrees C

## Choice of 3 wash programmes


- 1** Quick Cycle - only 35 minutes, including drying
- 2** Standard Cycle - 65 minutes, including drying
- 3** Test Cycle – for your weekly test regime, 50 minutes

## Stages within a wash cycle


1: Short Cycle

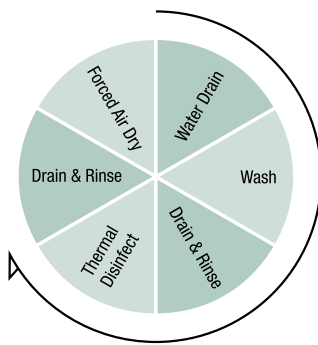




2: Standard Cycle

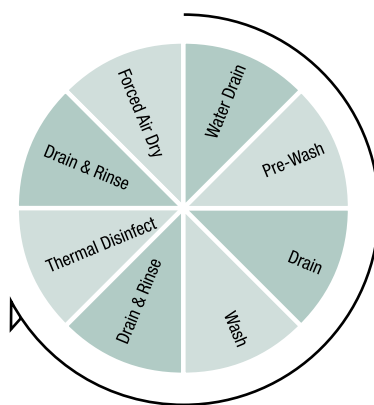

3: Test Cycle

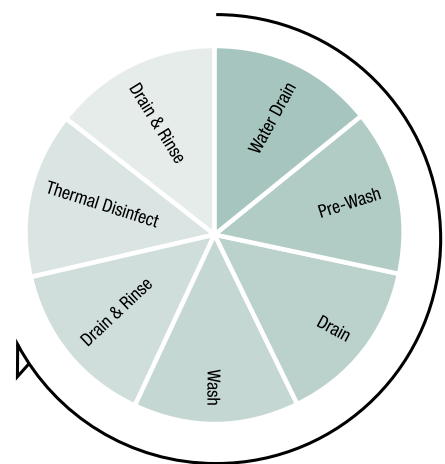

## Furniture and Accessories

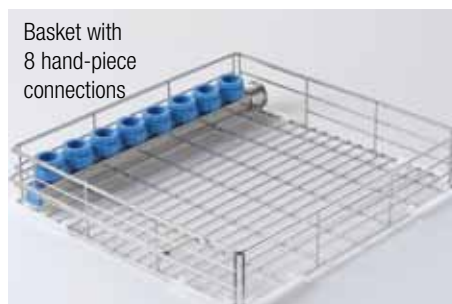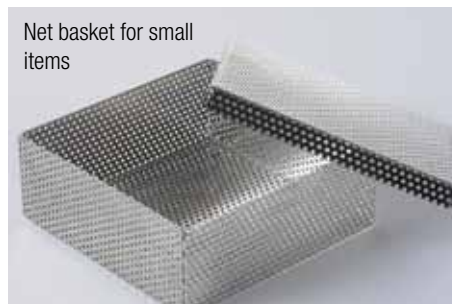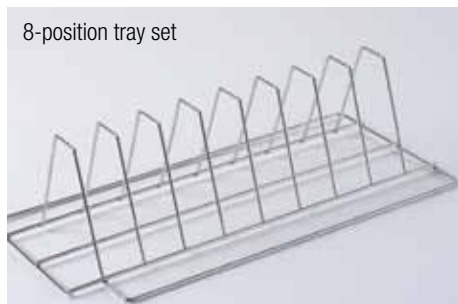

# Infection Control?

## 3 Easy steps to better practice with Prestige Medical

The UltraClean 3 bench top and under bench washer disinfectors are just one key step in your quest to comply with and exceed the guidelines for infection control.

Prestige Medical provides integrated decontamination solutions – products and services to help you at every step of your infection control procedure of washing, sterilizing and recording data.



## Sterilization as easy as 1-2-3

Enjoy the benefits of a fully integrated decontamination solution:

Prestige Medical products all work seamlessly together and are designed to give you a smooth, hands-free decontamination process.



**step 1**  
WASH  
Use the fully compliant UltraClean 3 washer disinfectant

**step 2**  
STERILIZE  
Use the fast effective Optima autoclave

**step 3**  
RECORD & TEST  
Conduct and record tests as specified by HTM 01-05

= full compliance  
To Best Practice Requirements



Contact us for advice on infection control, test procedures, decontamination rooms and equipment. [sales@prestigemedical.co.uk](mailto:sales@prestigemedical.co.uk) or telephone 01254 844 103.

# Prestige Medical - integrated decontamination solutions

For over 20 years, Prestige Medical has been respected as a class leader in the development and manufacture of sterilization and disinfection equipment. With over 1 million units delivered into over 80 countries worldwide, we have built an enviable reputation for innovation, performance and trust.

With infection control standards becoming increasingly stringent, a new breed of equipment is now required that is capable of meeting these challenges. The UltraClean 3 Washer Disinfector helps you surpass these standards. As part of our new generation of Integrated Decontamination Solutions we are able to take infection control to even greater levels. Our 3 step system has products and services that give you confidence at every stage of your decontamination and sterilization procedures.

Additionally, all Prestige Medical products and services are backed by a team of expert technicians, service engineers and customer support staff. Their mission is to deliver the most effective solution, impart their specialist knowledge of dental infection control and design a system that's exactly right for you.

## standards and guidelines

Designed and manufactured to comply with the following standards and guidelines:

- 93/42/EEC Medical Device Directive
- EN15883-1
- HTM 2030



## Contact us for advice on infection control



Prestige Medical Limited  
East House, Duttons Way,  
Shadsworth Business Park,  
Blackburn BB1 2QR

Tel : +44 (0) 1254 682 622  
Fax : +44 (0) 1254 682 606

[www.prestigemedical.co.uk](http://www.prestigemedical.co.uk)  
[sales@prestigemedical.co.uk](mailto:sales@prestigemedical.co.uk)

@PrestigeMedica1

Registered in England  
Reg No. 2826793  
Vat No. 633707835

Registered Office:  
Barlow Road,  
Aldermans Green Industrial Estate,  
Coventry CV2 2LD

Prestige Medical is a registered trademark



Due to a continuous process of development, Prestige Medical reserves the right to amend the specification of products and services without notice.