

Innovating Tradition.

About Future Lap[®]

The hybrid Future Lap[®] line consists of a high quality and easily detachable handle, which ensures comfortable and longterm repair free multiple use and a disposable tip-sheath assembly with all the benefits of single-use instruments.



Future Lap[®]



FUTURE LAP

How it



KEY FEATURES

- ✓ Hybrid instrument line with ball-socket-locking between handle and sheath
- ✓ Demountable into tip-sheath combination and handle
- ✓ New disposable instrument tip-sheath assembly for each procedure
 - ✓ Autoclavable handle part
 - ✓ Lengths of 330 and 450 mm on request
- ✓ No hidden spaces to accommodate blood and debris between reusable and disposable part
 - ✓ High-quality, ergonomic handles (System CS[®] or Endo Lock[®]) with 45° HF-port
- ✓ Quick, easy and safe assembly
- ✓ Combines benefits of single-use and reusable instruments
- ✓ Nearly no need for repair and maintenance



STEP

- ✓ Open the handle completely, until the socket coupling at the distal end becomes visible.



STEP

- ✓ Place the ball coupling of the tip into the sheath's inlet and close the handle completely to fix the tip-sheath combination.






STEP

- ✓ Turn the nut on the sheath counterwise onto the handle.





HYBRID LINE works

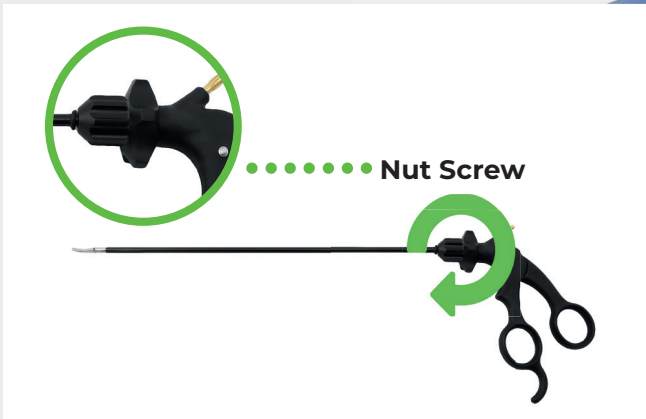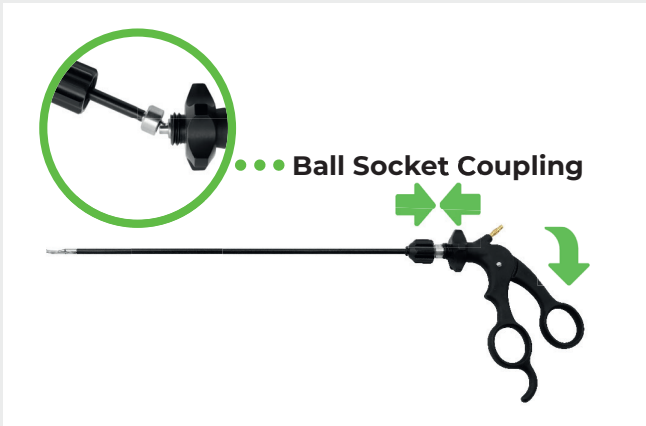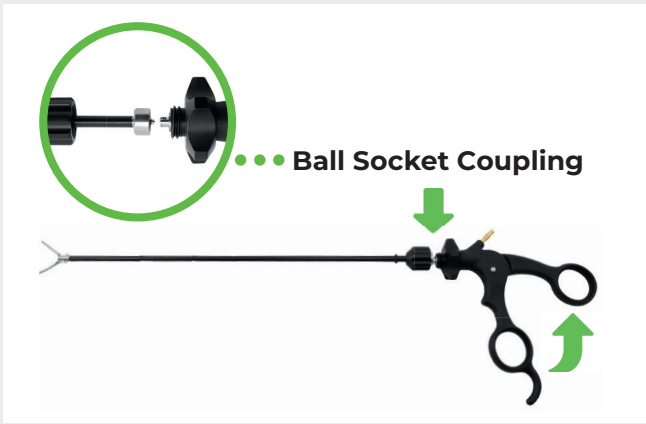



TOPSELLERS

Ordering Information

FUTURE LAP TIP-SHEATH SHEATH PART#

COMMON HANDLE

Fenestrated grasper
flat jaws
Ø 5 mm
length 330mm
UOM 10 pcs.



Polyamide handle
ergonomic
w/o ratchet
System CS[®]



Item # 13-1301DDI

Item # 13-1367CS

Maryland grasper
Ø 5 mm
length 330mm
UOM 10 pcs.



Polyamide handle
ergonomic
w/o ratchet
System CS[®]



Item # 13-1318DDI

Item # 13-1367CS

Endo clinch
Ø 5 mm
length 330mm
UOM 10 pcs.



Polyamide handle
ergonomic
universal ratchet
System CS[®]



Item # 13-1323DDI

Item # 13-1384CS

Metzenbaum scissors
curved, serrated
Ø 5 mm
length 330mm
UOM 10 pcs.



Polyamide handle
ergonomic
w/o ratchet
System CS[®]



Item # 13-1466DDI

Item # 13-1367CS

Metzenbaum scissors
Ø 5 mm
length 330mm
UOM 10 pcs.



Polyamide handle
ergonomic
w/o ratchet
System CS[®]



Item # 13-1468DDI

Item # 13-1367CS



www.ackermanninstrumente.de