




Suture Needles

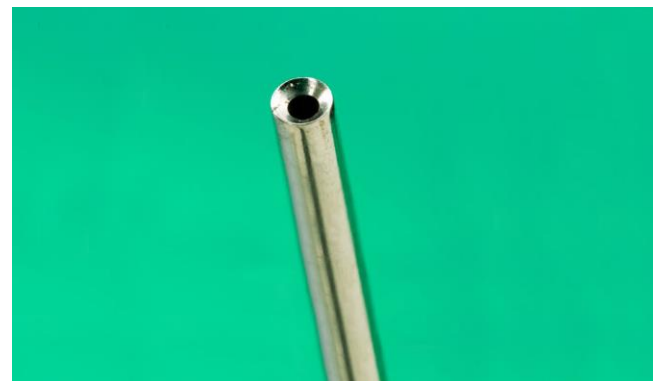

Barber of Sheffield, Unit 25 Shortwood Court, Shortwood Business Park,
Deame Valley Parkway, Barnsley, S74 9LH

Tel: +44 (0) 114 2576022 Fax: +44 (0) 114 2576222 Email: sales@barberofsheffield.co.uk

www.barberveterinary.co.uk

Barber of Sheffield Needles

Barber of Sheffield are the market leaders for eyed suture needles, with our needles supplied worldwide to veterinary practices, hospitals and clinics.

Manufactured from high grade stainless steel, our needles undergo a rigorous quality control procedure to ensure a high standard of sharpness, rigidity, durability and presentation. The needles are hand honed for added sharpness to ensure a smooth, less traumatic passage through tissue.

Eyed needles can be offered in regular, French or spring eye and in conventional cutting or reverse cutting. Conventional cutting needles have the cutting edge on the inside of the needle curvature and hence directed towards the wound, whereas reverse cutting needles have the cutting edge on the outside, with the flat surface directed towards the wound. Round bodied needles taper gradually to a point, whereas triangular bodies have cutting edges along three sides.


We have a variety of suture needle shapes available, including 3/8 circle (curved), 1/2 curved, 1/2 circle and straight, as well as different body types - round and triangular. In addition to the standard suture needles, we also offer a post mortem needle range, which are larger and thicker. These are available in 3/8 circle (curved), 1/2 curved, straight and serpentine (s-shaped).

Barber can offer a range of special needles, such as drilled end and micro needles. For information on bespoke needles please contact us at sales@barberneedles.com

For further information on our needles, please contact the Medical & Veterinary Sales Manager:

Name: Ian Couldwell
Telephone: +44 (0) 114 3493797
Email: ian@barberofsheffield.co.uk

You can also find further information on all the Barber of Sheffield products on our 3 websites:

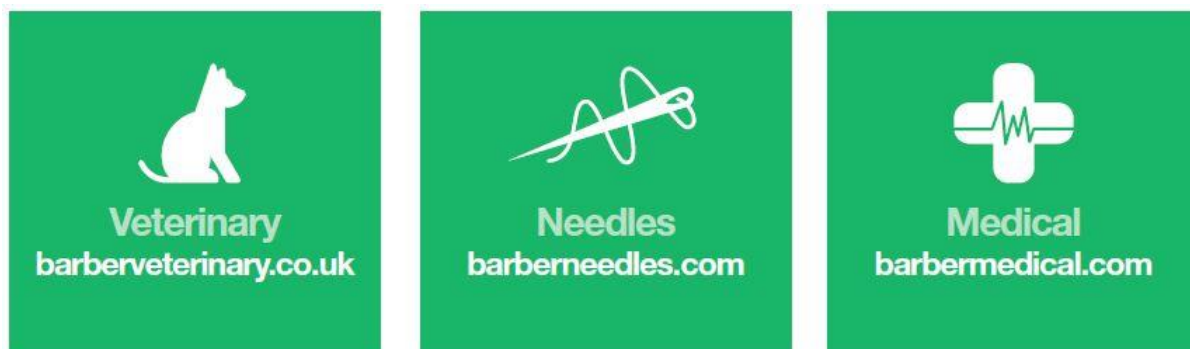

Barber of Sheffield, Unit 25 Shortwood Court, Shortwood Business Park,
Dearne Valley Parkway, Barnsley, S74 9LH

Tel: +44 (0) 114 2576022 Fax: +44 (0) 114 2576222 Email: sales@barberofsheffield.co.uk

www.barbveterinary.co.uk

Needle Eye Types

Regular Eye



Spring Eye



Needle Shapes




Straight



Half Curved



3/8 Circle

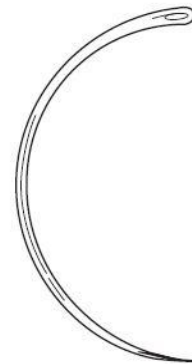

1/2 Circle

Needle Point Types



Triangular Cutting
Point



Round
Point



Trocar
Point

Barber of Sheffield, Unit 25 Shortwood Court, Shortwood Business Park,
Dearne Valley Parkway, Barnsley, S74 9LH

Tel: +44 (0) 114 2576022 Fax: +44 (0) 114 2576222 Email: sales@barberofsheffield.co.uk


www.barberveterinary.co.uk

3/8 Circle Eyed Suture Needles

The 3/8 circle needles, which are also referred to as 'curved', are the most commonly used needle. The curvature makes it easy to manipulate in large or superficial wounds. The needles are available in packs of 12, in the following variants:

Product Code	Eye Type?	Point Type?	Body Type
NCTC.01 to NCTC.20	Regular	Triangular Cutting	Conventional
NRBC.01 to NRBC.20	Regular	Round Taper	Conventional
REVNCTC.01 to REVNCTC.20	Regular	Triangular Cutting	Reverse
SCTC.01 to SCTC.20	Spring Eye	Triangular Cutting	Conventional
SRBC.01 to SRBC.20	Spring Eye	Round Taper	Conventional
REVSCTC.01 to REVSCTC.20	Spring Eye	Triangular Cutting	Reverse

The needles are available in sizes 1 to 20 (80mm to 23mm), which are detailed on the diagram below:

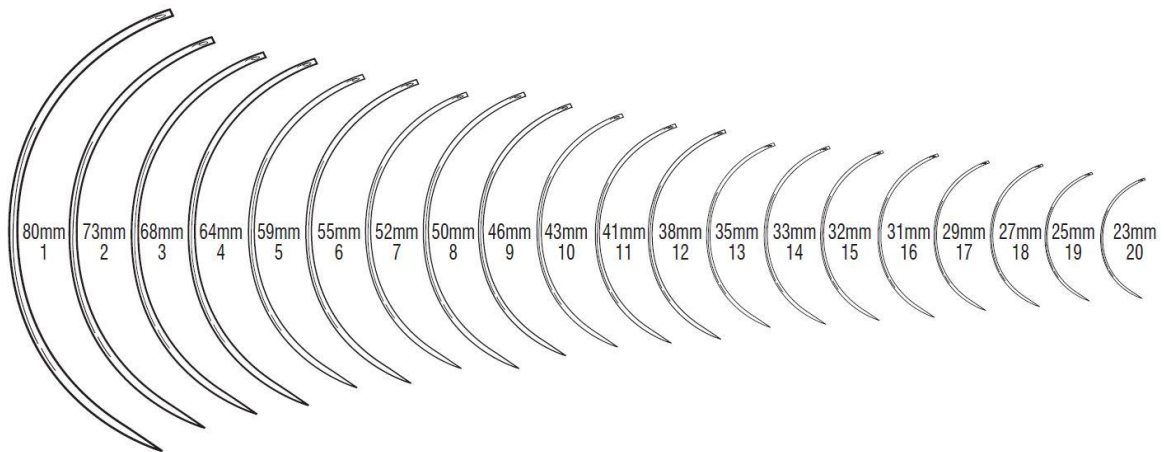

Barber of Sheffield, Unit 25 Shortwood Court, Shortwood Business Park,
Dearne Valley Parkway, Barnsley, S74 9LH

Tel: +44 (0) 114 2576022 Fax: +44 (0) 114 2576222 Email: sales@barberofsheffield.co.uk


www.barberveterinary.co.uk

1/2 Curved Eyed Suture Needles

The 1/2 curved needles are commonly used for the laparoscopic surgical technique and are available in packs of 12, in the following variants:

Product Code	Eye Type?	Point Type?	Body Type
NHCTC.01 to NHCTC.20	Regular	Triangular Cutting	Conventional
NHRB.01 to NHRB.20	Regular	Round Taper	Conventional
SHCTC.01 to SHCTC.20	Spring Eye	Triangular Cutting	Conventional
SHRB.01 to SHRB.20	Spring Eye	Round Taper	Conventional

The needles are available in sizes 1 to 20 (80mm to 23mm), which are detailed on the diagram below:

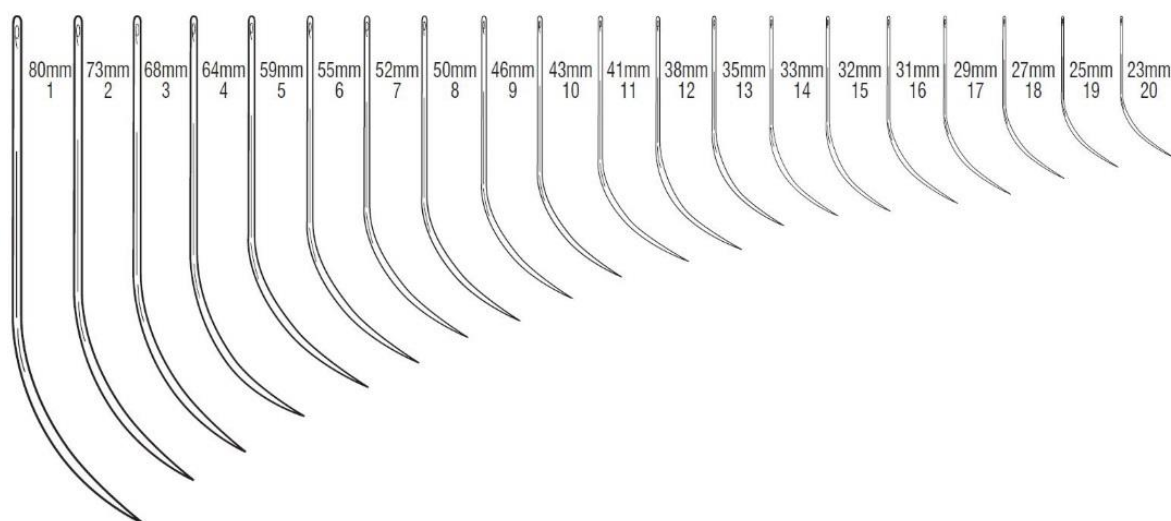

Barber of Sheffield, Unit 25 Shortwood Court, Shortwood Business Park,
Dearne Valley Parkway, Barnsley, S74 9LH

Tel: +44 (0) 114 2576022 Fax: +44 (0) 114 2576222 Email: sales@barberofsheffield.co.uk


www.barberveterinary.co.uk

1/2 Circle Eyed Suture Needles

The 1/2 circle needles are usually used for surgery in confined locations and are available in packs of 12, in the following variants:

Product Code	Eye Type?	Point Type?	Body Type
NHCRTC.01 to NHCRTC.20	Regular	Triangular Cutting	Conventional
NHCRRB.01 to NHCRRB.20	Regular	Round Taper	Conventional
REVNHCRTC.01 to REVNHCRTC.20	Regular	Triangular Cutting	Reverse
SHCRTC.01 to SHCRTC.20	Spring Eye	Triangular Cutting	Conventional
SHCRRB.01 to SHCRRB.20	Spring Eye	Round Taper	Conventional

The needles are available in sizes 1 to 20 (80mm to 23mm), which are detailed on the diagram below:

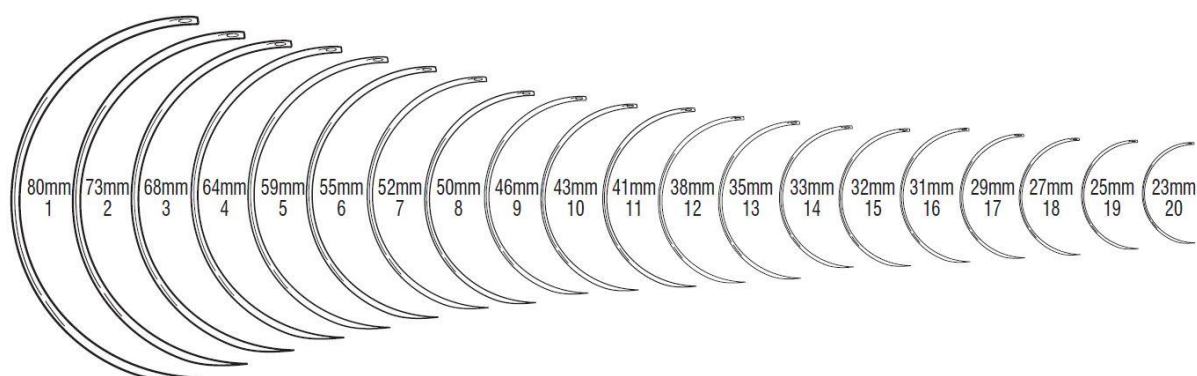

Barber of Sheffield, Unit 25 Shortwood Court, Shortwood Business Park,
Dearne Valley Parkway, Barnsley, S74 9LH

Tel: +44 (0) 114 2576022 Fax: +44 (0) 114 2576222 Email: sales@barberofsheffield.co.uk


www.barberveterinary.co.uk

Straight Eyed Suture Needles

The straight needles are used when the tissue is easily accessible, so they are usually used for skin closure or nerve/vessel repair. The needles are available in packs of 12, in the following variants:

Product Code	Eye Type?	Point Type?	Body Type
NSTC.01 to NSTC.20	Regular	Triangular Cutting	Conventional
NSRB.01 to NSRB.20	Regular	Round Taper	Conventional

The needles are available in sizes 1 to 20 (80mm to 23mm), which are detailed on the diagram below:

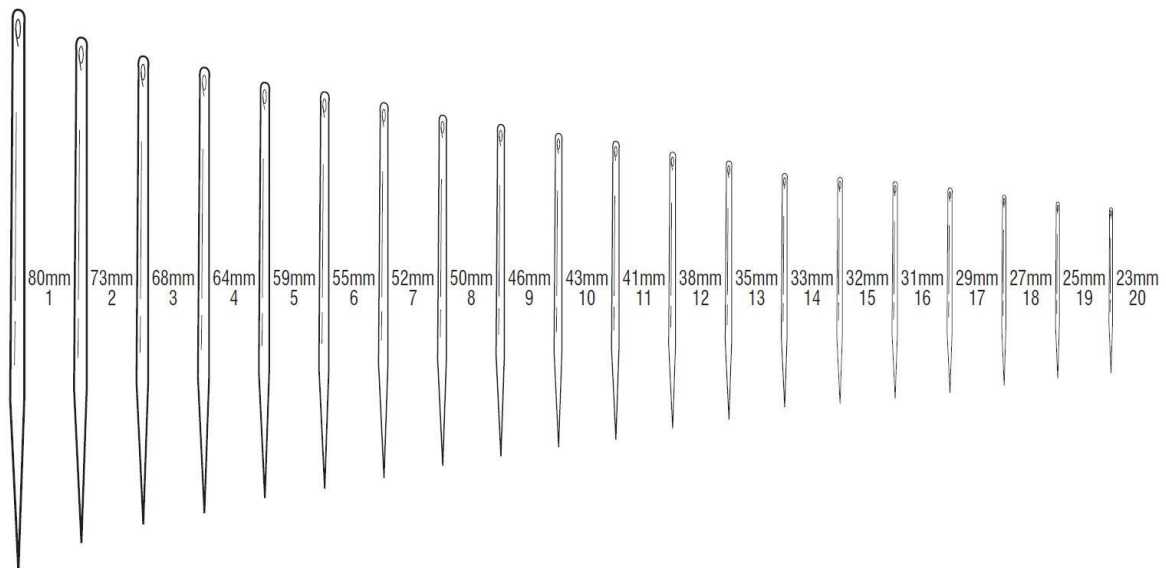

Barber of Sheffield, Unit 25 Shortwood Court, Shortwood Business Park,
Dearne Valley Parkway, Barnsley, S74 9LH

Tel: +44 (0) 114 2576022 Fax: +44 (0) 114 2576222 Email: sales@barberofsheffield.co.uk

www.barberveterinary.co.uk

Mayo 1/2 Circle Catgut Taper Needles


The Mayo 1/2 circle catgut taper needles are usually used for tissue closure surgery. The needles can be supplied non-sterile or sterile.

Taper Point – Non Sterile

The following non-sterile taper point needles are supplied in packs of 12:

Product Code	Eye Type?	Point Type?	Body Type
NMHCTTT.01 to NMHCTTT.06	Regular	Taper	Conventional

The needles are available in sizes 1 to 6 (51mm to 25mm), which are detailed on the diagram below:

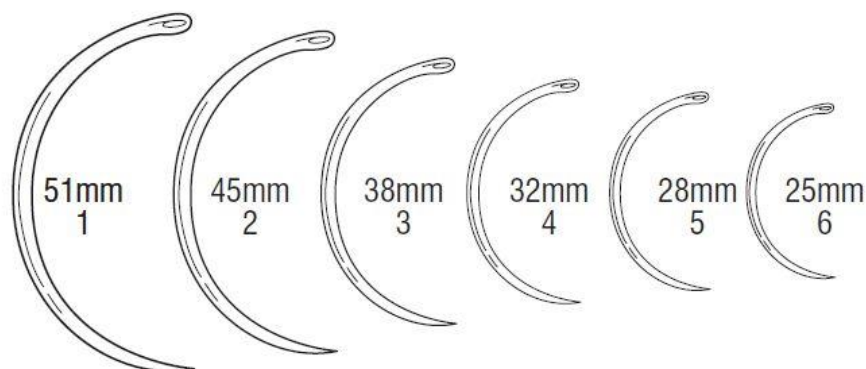

Taper Point – Sterile

The sterile taper point needles are a slightly different size and are supplied in sterile pouches of 2 needles per pouch. Each of the below codes is sold as a box of 25 sterile pouches:

Product Code	Size	Eye Type?	Point Type?	Body Type
824-S-2	45mm	Regular	Taper	Conventional
824-S-3	38mm	Regular	Taper	Conventional
824-S-4	35mm	Regular	Taper	Conventional
824-S-5	32mm	Regular	Taper	Conventional
824-S-6	28mm	Regular	Taper	Conventional

Barber of Sheffield, Unit 25 Shortwood Court, Shortwood Business Park,
Dearne Valley Parkway, Barnsley, S74 9LH

Tel: +44 (0) 114 2576022 Fax: +44 (0) 114 2576222 Email: sales@barberofsheffield.co.uk

www.barberveterinary.co.uk

Mayo 1/2 Circle Catgut Trocar Needles


The Mayo 1/2 circle catgut trocar needles are usually used for tissue closure surgery. The needles can be supplied non-sterile or sterile.

Trocar Point – Non Sterile

The following non-sterile trocar point needles are supplied in packs of 12:

Product Code	Eye Type?	Point Type?	Body Type
NMHCR.T.01 to NMHCR.T.06	Regular	Trocar	Conventional

The needles are available in sizes 1 to 6 (51mm to 25mm), which are detailed on the diagram below:

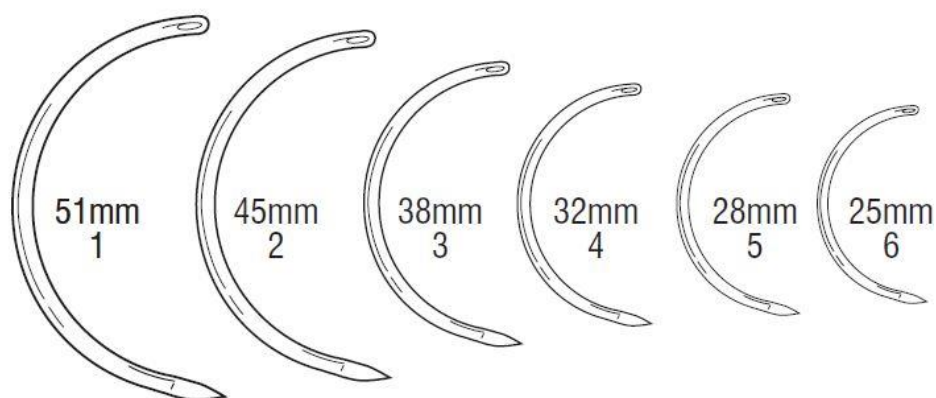

Trocar Point – Sterile

The sterile trocar point needles are a slightly different size and are supplied in sterile pouches of 2 needles per pouch. Each of the below codes is sold as a box of 25 sterile pouches:

Product Code	Size	Eye Type?	Point Type?	Body Type
826-S-1	51mm	Regular	Trocar	Conventional
826-S-2	45mm	Regular	Trocar	Conventional
826-S-3	38mm	Regular	Trocar	Conventional
826-S-4	35mm	Regular	Trocar	Conventional
826-S-5	32mm	Regular	Trocar	Conventional
826-S-6	28mm	Regular	Trocar	Conventional

Barber of Sheffield, Unit 25 Shortwood Court, Shortwood Business Park,
Dearne Valley Parkway, Barnsley, S74 9LH

Tel: +44 (0) 114 2576022 Fax: +44 (0) 114 2576222 Email: sales@barberofsheffield.co.uk

www.barberveterinary.co.uk

Post Mortem Needles

Our Barber Medical department is the market leading mortuary/autopsy consumables supplier in the UK. Our range of post mortem needles have been specifically designed for the needs of mortuaries, but due to their extra thickness and larger sizes, these are also used across various industries, including veterinary. The post mortem needles are available in 3/8 circle (curved), 1/2 curved, serpentine (s-shape), 1/2 Circle and straight, in packs of 12 and in the following variants:

3/8 Circle

Product Code	Needle Length	Eye Type?	Point Type?
NPOSTC03	3inch	Regular	Triangular Cutting
NPOSTC04	4inch	Regular	Triangular Cutting
NPOSTC05	5inch	Regular	Triangular Cutting
NPOSTC06	6inch	Regular	Triangular Cutting

1/2 Curved

Product Code	Needle Length	Eye Type?	Point Type?
NPOSTHCR.03	3inch	Regular	Triangular Cutting
NPOSTHCR3.5INCH	3.5inch	Regular	Triangular Cutting
NPOSTHCR04	4inch	Regular	Triangular Cutting
NPOSTHCR045	4.5inch	Regular	Triangular Cutting
NPOSTHCR05	5inch	Regular	Triangular Cutting
NPOSTHCR055	5.5inch	Regular	Triangular Cutting
NPOSTHCR.06	6inch	Regular	Triangular Cutting

Serpentine

Product Code	Needle Length	Eye Type?	Point Type?
NPOSTSER5	2.25inch	Regular	Triangular Cutting
NPOSTSER4	2.5inch	Regular	Triangular Cutting
NPOSTSER3	3.25inch	Regular	Triangular Cutting
NPOSTSER2	4inch	Regular	Triangular Cutting
NPOSTSER01	5inch	Regular	Triangular Cutting
NPOSTSER0	5.5inch	Regular	Triangular Cutting
NPOSTSER6INCH	6inch	Regular	Triangular Cutting

1/2 Circle

Product Code	Needle Length	Eye Type?	Point Type?
NPOSTHCIR4	4inch	Regular	Triangular Cutting
NPOSTHCR3/0	4.5inch	Regular	Triangular Cutting

Barber of Sheffield, Unit 25 Shortwood Court, Shortwood Business Park,
Dearne Valley Parkway, Barnsley, S74 9LH


Tel: +44 (0) 114 2576022 Fax: +44 (0) 114 2576222 Email: sales@barberofsheffield.co.uk

www.barberveterinary.co.uk

Straight

Product Code	Needle Length	Eye Type?	Point Type?
NPOSTSTR.04	4inch	Regular	Triangular Cutting
NPOSTSTR05	5inch	Regular	Triangular Cutting
NPOSTSTR.06	6inch	Regular	Triangular Cutting
NPOSTSTR.09	9inch	Regular	Triangular Cutting

The majority of the post mortem needle sizes and shapes available are detailed on the diagram below:

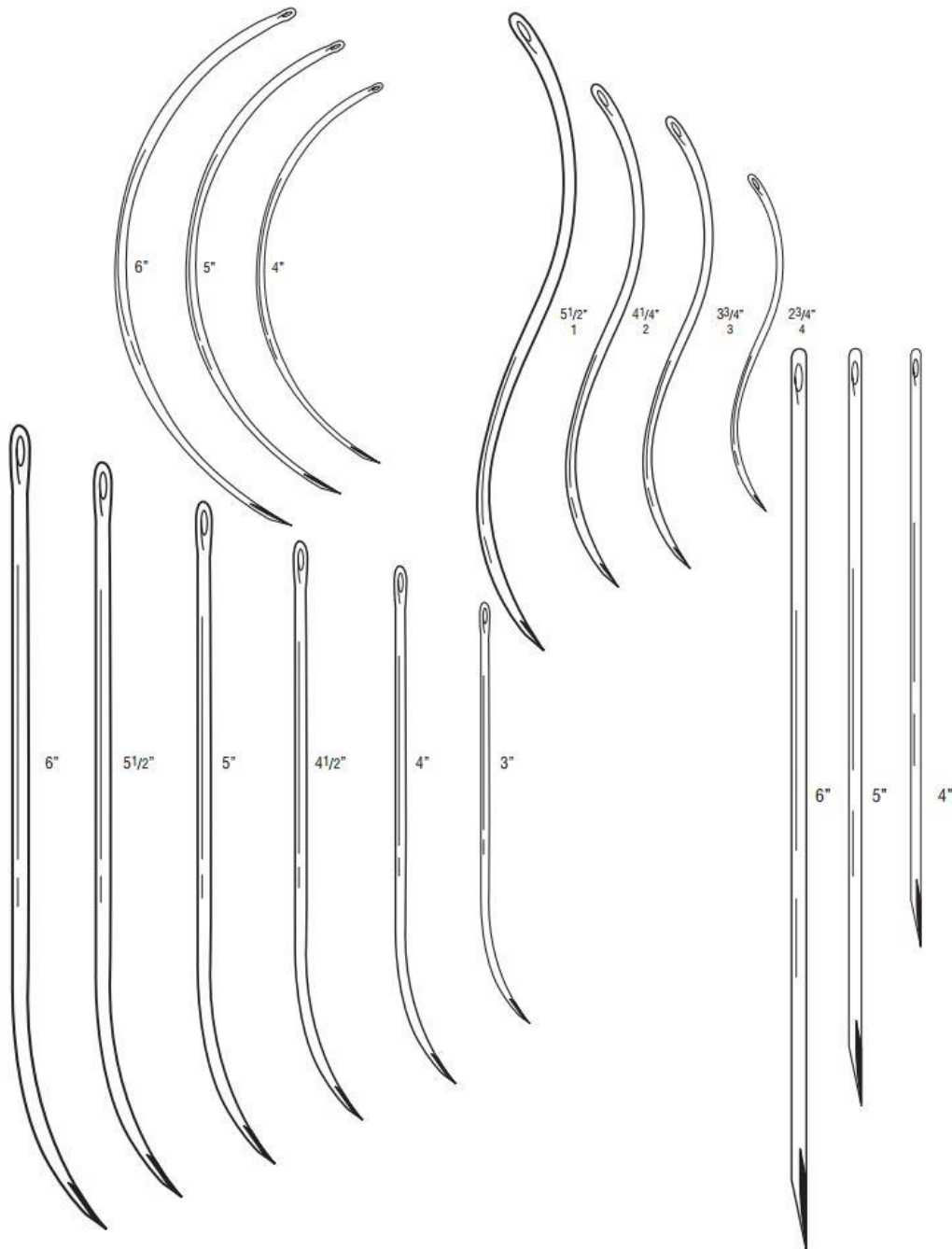

Barber of Sheffield, Unit 25 Shortwood Court, Shortwood Business Park,
Dearne Valley Parkway, Barnsley, S74 9LH

Tel: +44 (0) 114 2576022 Fax: +44 (0) 114 2576222 Email: sales@barberofsheffield.co.uk


www.barberveterinary.co.uk

Loopuyt Needles

The Loopuyt needles are often used in by mortuaries for post mortems, but can be used for other applications. These needles have a double curve for ease of grip and are sold in packs of 12:

Product Code	Needle Length	Eye Type?	Point Type?
LOOPUYTS1	80mm	Regular	Triangular Cutting
LOOPUYTS2	72mm	Regular	Triangular Cutting
LOOPUYTS68MM	68mm	Regular	Triangular Cutting
LOOPUYTS60MM	60mm	Regular	Triangular Cutting
LOOPUYTS3	55mm	Regular	Triangular Cutting
LOOPUYTS4	52mm	Regular	Triangular Cutting
LOOPUYTS5	46mm	Regular	Triangular Cutting

The Loopuyt needles are standardly available in sizes 1-5, but we have added 2 additional sizes (60mm & 68mm):

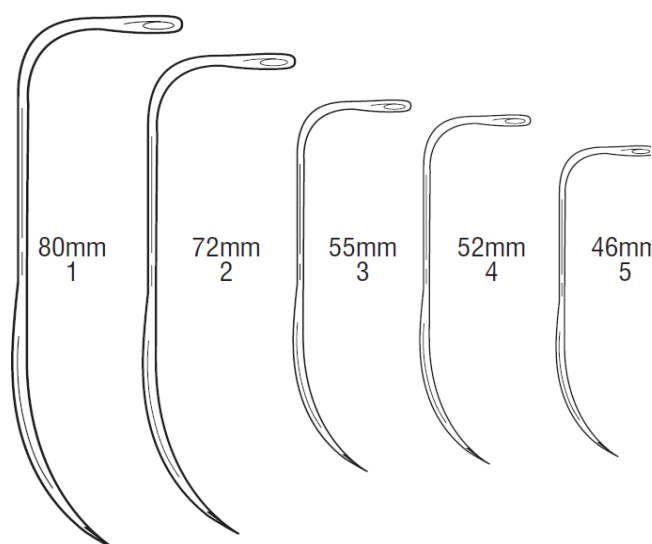

Colt's Tension Needles

The colt's tension needles are 3/8 circle (curved) and are used for wound closure. These needles are available in 2 sizes:

Product Code	Needle Length	Eye Type?	Point Type?
NCOLTS3	3inch	Regular	Triangular Cutting
NCOLTS4	4inch	Regular	Triangular Cutting

Barber of Sheffield, Unit 25 Shortwood Court, Shortwood Business Park,
Deame Valley Parkway, Barnsley, S74 9LH


Tel: +44 (0) 114 2576022 Fax: +44 (0) 114 2576222 Email: sales@barberofsheffield.co.uk

www.barberveterinary.co.uk

Keith's Abdominal Needles

The Keith abdominal needles are straight needles, which are used for suturing abdominal tissue during veterinary procedures and are available in packs of 12, in the following sizes:

Product Code	Needle Length	Eye Type?	Point Type?
KEITH11/2INCHDOZ824-11/2	1.5inch	Regular	Triangular Cutting
KEITH1.75INCH824-13/4DOZ	1.75inch	Regular	Triangular Cutting
KEITH2INCHDOZ	2inch	Regular	Triangular Cutting
KEITH21/2INCH824-21/2DOZ	2.5inch	Regular	Triangular Cutting
KEITH2.75INCHDOZ	2.75inch	Regular	Triangular Cutting
KEITH3INCHDOZ	3inch	Regular	Triangular Cutting
KEITH33/4INCHDOZ	3.75inch	Regular	Triangular Cutting

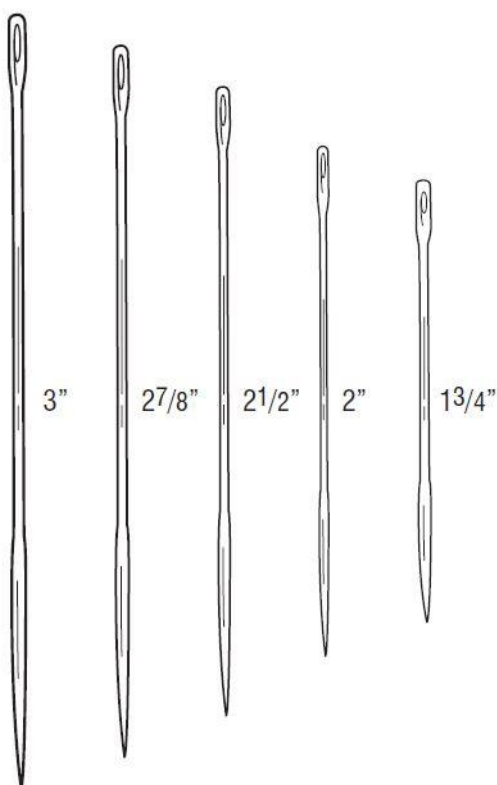

Barber of Sheffield, Unit 25 Shortwood Court, Shortwood Business Park,
Dearne Valley Parkway, Barnsley, S74 9LH

Tel: +44 (0) 114 2576022 Fax: +44 (0) 114 2576222 Email: sales@barberofsheffield.co.uk

www.barberveterinary.co.uk

Ferguson Needles

The Ferguson needles are 1/2 circle taper point needles, which are commonly used for suturing organs and vessels in veterinary procedures and are available in packs of 12, in the following sizes:

Product Code	Eye Type?	Point Type?	Body Type?
NFERGUSON.03	Regular	Round Taper	Conventional
NFERGUSON.04	Regular	Round Taper	Conventional
NFERGUSON.06	Regular	Round Taper	Conventional
NFERGUSON.08	Regular	Round Taper	Conventional
NFERGUSON.09	Regular	Round Taper	Conventional
NFERGUSON.10	Regular	Round Taper	Conventional
NFERGUSON.11	Regular	Round Taper	Conventional
NFERGUSON.12	Regular	Round Taper	Conventional
NFERGUSON.13	Regular	Round Taper	Conventional

Cruciate Needles


The cruciate needles are used during veterinary cruciate repair procedures. We offer these in both non-sterile and sterile with leaderline options. The needles are sold in packs of 6 and are available in the following variants:

Non-Sterile

Product Code	Cruciate Size	Eye Type?	Point/Body Type?
CRUCIATEVERYSMALL	Very Small	Regular	Round Taper
CRUCIATESMALL	Small	Regular	Round Taper
CRUCIATEMED	Medium	Regular	Round Taper
CRUCIATELARGE	Large	Regular	Round Taper

Sterile

Product Code	Cruciate Size	Eye Type?	Point/Body Type?
CRUCIATE5STERILE	Very Small	Regular	Round Taper
CRUCIATE4STERILE	Small	Regular	Round Taper
CRUCIATE3STERILE	Medium	Regular	Round Taper

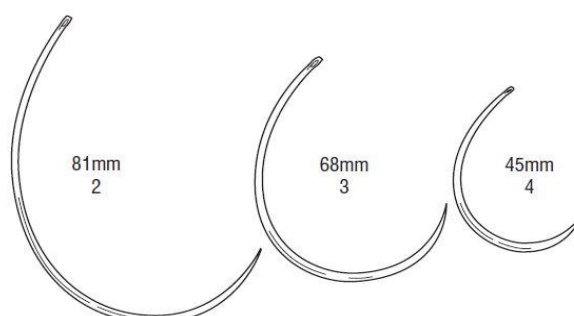

Barber of Sheffield, Unit 25 Shortwood Court, Shortwood Business Park,
Dearne Valley Parkway, Barnsley, S74 9LH

Tel: +44 (0) 114 2576022 Fax: +44 (0) 114 2576222 Email: sales@barberofsheffield.co.uk

www.barberveterinary.co.uk

Hypodermic Needles



Hypodermic disposable needles are a hollow needle commonly used with a syringe to inject substances into the body or extract fluids from it.

We sell both the sterile hypodermic single use needles and the hypodermic Luer & Record re-usable needles:

Sterile - Single Use

The single use sterile hypodermic needles are sold in packs of 100 in the following sizes:

Product Code	Gauge	Length (Inch)
NEEDLEDISP14GX1.5INCH	14G	1.5
NEEDLEDISP16GX1.5INCH	16G	1.5
NEEDLEDISP16GX1INCH	16G	1
NEEDLEDISP17GX0.5INCH	17G	0.5
NEEDLEDISP17GX0.75	17G	0.75
NEEDLEDISP17GX1INCH	17G	1
NEEDLEDISP18GX0.5INCH	18G	0.5
NEEDLEDISP18GX0.75INCH	18G	0.75
NEEDLEDISP18GX1INCH	18G	1
NEEDLEDISP18GX1.5INCH	18G	1.5
NEEDLEDISP19GX1INCH	19G	1
NEEDLEDISP20GX0.5INCH	20G	0.5
NEEDLEDISP21GX5/8INCH	21G	0.62

Barber of Sheffield, Unit 25 Shortwood Court, Shortwood Business Park,
Dearne Valley Parkway, Barnsley, S74 9LH

Tel: +44 (0) 114 2576022 Fax: +44 (0) 114 2576222 Email: sales@barberofsheffield.co.uk

www.barberveterinary.co.uk

Luer & Record Needles

The re-usable stainless steel needles are available with a Luer or Record hub. The re-usable needles are either sold as singles or packs of 12

Luer Hub:

Product Code	Gauge	Length (Inch)	Pack Size
LUER13GX1.75SINGLE"	13G	1.75	1
LUER13GX212'S"	13G	2	12
LUER14GX5/812'S"	14G	0.62	12
LUER14GX3/412"	14G	0.75	12
LUER14GX1.2512'S"	14G	1.25	12
LUER14GX1.512'S"	14G	1.5	12
LUER14GX4INCH	14G	4	1
LUER16GX3/812"	16G	0.37	12
LUER16GX1/212"	16G	0.5	12
LUER16GX3/412"	16G	0.75	12
LUER16GX112"	16G	1	12
LUER16GX1.2512"	16G	1.25	12
LUER16GX1.512"	16G	1.5	12
LUER17GX1/212"	17G	0.5	12
LUER17GX3/412"	17G	0.75	12
LUER17GX112'S"	17G	1	12
LUER17GX1.512'S	17G	1.5	12
LUER17GX3INCH	17G	3	1
LUER17GX4.07INCH	17G	4.07	1
LUER18GX1/212"	18G	0.5	12
LUER18GX3/412"	18G	0.75	12
LUER18GX112"	18G	1	12
LUER18GX1.2512S"	18G	1.25	12
LUER19GX3/812'S"	19G	0.37	12
LUER19GX1/212"	19G	0.5	12
LUER19GX3/412"	19G	0.75	12
LUER19GX112'S"	19G	1	12
LUER19GX3INCH	19G	3	1
LUER19GX4.07INCH	19G	4.07	1
LUER20GX1/212"	20G	0.5	12
LUER20GX5/812'S"	20G	0.62	12
LUER20GX3/412'S"	20G	0.75	12
LUER20G1.512"	20G	1.5	12
LUER21GX3/812"	21G	0.37	12
LUER21GX1/212"	21G	0.5	12
LUER22GX1/412"	22G	0.25	12
LUER22GX3/812'S"	22G	0.37	12
LUER22GX1/212'S	22G	0.5	12
LUER23GX212'S"	23G	2	12

Barber of Sheffield, Unit 25 Shortwood Court, Shortwood Business Park,
Dearne Valley Parkway, Barnsley, S74 9LH

Tel: +44 (0) 114 2576022 Fax: +44 (0) 114 2576222 Email: sales@barberofsheffield.co.uk

www.barberveterinary.co.uk

Record Hub:

Product Code	Gauge	Length (Inch)	Pack Size
RECORD13GX13/4"	13G	1.75	1
RECORD13GX212'S"	13G	2	12
RECORD14GX1.512'S"	14G	1.5	12
RECORD14GX1.7512"	14G	1.75	12
RECORD14GX212'S"	14G	2	12
RECORD15GX112'S"	15G	1	12
RECORD16GX1/212'S"	16G	0.5	12
RECORD16GX1.512'S"	16G	1.5	12
RECORD18GX1/212'S"	18G	0.5	12
RECORD19GX1/212'S"	19G	0.5	12
RECORD19GX3/412'S"	19G	0.75	12
RECORD20GX3/812'S	20G	0.37	12
RECORD20GX1/212'S"	20G	0.5	12

Barber of Sheffield, Unit 25 Shortwood Court, Shortwood Business Park,
Dearne Valley Parkway, Barnsley, S74 9LH

Tel: +44 (0) 114 2576022 Fax: +44 (0) 114 2576222 Email: sales@barberofsheffield.co.uk

www.barberveterinary.co.uk